

845 RANGE OF MOBILE COMMODES

Model Types

Chrome plated steel adjustable height mobile commode chairs, c/w swinging drop down & detachable armrests, detachable back, 4 x 5" brake castors, plastic commode seat with square bowl, upholstered overseat & backrest, built in pushing handle.

84518	18" seat width	Weight capacity	190 Kgs (30 st)
84520	20" seat width	Weight capacity	190 Kgs (30 st)
84522	22" seat width	Weight capacity	190 Kgs (30 st)



Operation

- ☛ **WARNING!** - These products should always be used under supervision.

Brakes/Castors

The Days Healthcare 845 range of mobile commodes is fitted with 4 fully braked castors. Each castor can be braked by pushing down the lever on the castor. This prevents the wheels from rotating or swivelling. To release the brake the lever should be pushed upward.

- ☛ **WARNING!** – Brakes should always be applied with the castors facing forward (fig. 1) when entering or leaving the chair. With the brakes applied the commode may still slide on a smooth floor or if it is tipped forward by the user. Supervision is essential while getting on or off the commode.

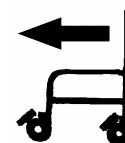


fig. 1


Bowl Removal

The commode bowl is removed by sliding out from the back of the seat.

Armrests

The 845 models have two armrests, which slot into brackets by the side of the chair. The arms can be removed by pressing the lever towards the arm, raising the front side of the arm and then releasing the arm at the rear by pressing on the spring button. When re-fitting the arm, ensure that the spring button is locked in place. To lower the armrests only, press the lever towards the arm, raise the arm and then lower to one side until it is resting on the support bar.

- ☛ **WARNING!**
 - Ensure the arms are securely slotted into the brackets before use.
 - Do not attempt to lift the commode by the arms.

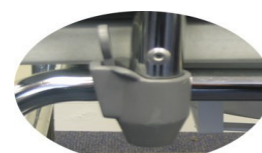

845PL002
Armrest Release Clip

Footrests

The 845 range of mobile commodes are fitted with hook on footrests as standard (Days Healthcare Part No. 845SA006). To adjust height of footrests, loosen the hexagonal bolt on bottom of tube until the footrest can move smoothly, slide down until the height required and tighten the hexagonal bolt until the footrest is rigid. For access, the footrest plates can be swivelled up vertically.



845SA006 footrests

Fitting of Detachable Backrest

The backrest is fitted by inserting the back into the frame and tightening up the screw with the hand knob, ensuring that the screw head is recessed in the hexagonal shape on the inside of the tube and the hand knob and washer is fitted on the outside.




845HB001 Hand Knob

Perfection Pan & Pan Rack (optional extra)

All commodes can be fitted with Perfection Pan Rack (part number ACC/2).

To fit the pan rack slide it onto the runners under the seat (SA154) and hook the rear crossbar over the plastic clip at the back of the commode seat. The rack will take the Perfection Pan (part number 640) and a disposable pan.

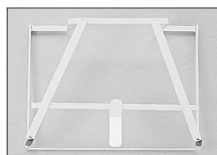

ACC/2 Perfection Pan Rack



SA154 Seat Assembly



No. 640 Bed Pan

Height Adjustment

To adjust the height take the load off the castors then press in the two buttons that lock-in the extension legs. Turn the leg slightly then pull the leg to the new leg setting. The two buttons should pop out fully to secure the leg. Make sure all four legs are set at the same height.

Maintenance

The Days Healthcare 845 Mobile Commodes are designed for minimum maintenance and easy, modular replacement of parts for long life. It is manufactured from chrome plated steel.

To clean, wipe over with a soft cloth. Use general purpose neutral detergent and warm water.

Do not use harsh abrasives.

To disinfect clean as above and wipe over with solution of one thousand parts per million available chlorine disinfectant (haztabs/precept), rinse and dry. Do not use any decontamination process that involves immersion.

Daily, before use, check the effective operation of the brakes and the secure fixing of the castors to the legs. Check that the commode is stable, i.e. all four castors are touching the floor and none of the legs are bent. Check that the commode arms are located correctly.

The chair is designed to allow damaged or broken parts to be changed easily. If any parts are damaged please contact the supplier of the chair or Days Healthcare. direct for replacement parts and fitting instructions.

Spare Parts

Square Pot	PL512
Blue Upholstered Overseat	18" (84518SA003), 20" (84520SA003), 22" (84522SA003)
Blue Upholstered Backrest	18" (84518SA004), 20" (84520SA004), 22" (84522SA004)
Brake Castors	845SA005
Spring Buttons to adj. ht.	SA001
Arms c/w fittings	845SA008(R) or 845SA009(L)
Blue Armpad only	845PU001

Warranty Information

The Days Healthcare 845 Mobile Commodes are warranted for one year from the date of purchase against faulty workmanship or materials. Please contact your supplier/dealer should a fault occur.

The warranty does not extend to the consequential costs resulting from fault clearance, in particular freight and travel costs, loss of earnings, expenses, etc.

The manufacturer will not accept responsibility for any damage or injury caused by misuse or non-observance of the instructions set out above.