

## SMART-RAIL

Sitting up and transferring out of bed is made easier with the innovative Smart-Rail. Unlike other bed rails, the Smart-Rail has two positions: closed and open. In the closed position (parallel to the side of the bed), the rail provides support for making adjustments while in bed. In the open position, the rail pivots out from the bed and touches down to the ground to assist movement to a standing position. The Smart-Rail is also easily installed between your mattress and box spring.




### SPECS

Rail (grip) Diameter	1.25"
Rail Width	16.5"
Rail Height	31.5"
Material	Alloy steel
Finish	Anti-microbial powder coat white (PCW)
Weight Capacity	300lbs
Code	SR-S <i>Smart-Rail</i>

Install  
Yourself

### PRODUCT DIMENSIONS

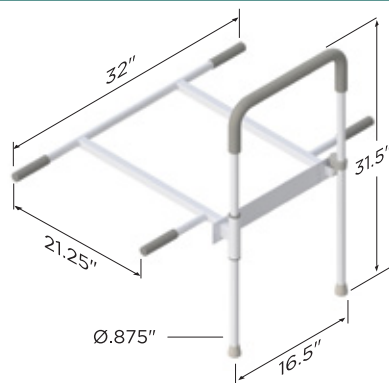

### SMART RAIL INSTALLATION ACCESSORIES

**Smart Rail Leg Extender** Increases standard rail height of 31.5" up to 34.5" in the open position

Code SR-LE *Smart-Rail Leg Extenders*

