

## Cushioning Products

# Standard Series

The original ROHO cushions ushered in a new era of skin protection and injury prevention, and today they still have no equal.

The DRY FLOATATION® technology redistributes pressure while decreasing friction and shear. Every ROHO Standard Series cushion is fully adjustable to fit your shape, minimize tissue deformation and facilitate blood flow.

ROHO HIGH PROFILE® Single Compartment Cushion

ROHO MID PROFILE™ Single Compartment Cushion

ROHO LOW PROFILE® Single Compartment Cushion



Pressure injury prevention  
and treatment



Adjustable



Lightweight



Washable



## ROHO® HIGH PROFILE® Single Compartment Cushion

The tall, 4.25 in. (10.5 cm) individual cells envelop you in air to redistribute pressure and decrease friction and shear. With a simple inflation adjustment, your ROHO will fit your exact shape to minimize tissue deformation and maximize blood flow.



Pressure injury prevention and treatment



Lightweight



Positioning



Washable



Adjustable

### Dimensions

Width: 11.5 in. to 26.75 in. / 29 cm to 68 cm

Depth: 11.5 in to 20 in. / 29 cm to 50.5 cm

Height: 4.25 in. / 10.5 cm

Weight: Approx. 3.25 lbs. / 1.5 kg

### Weight limit

No weight limit. However, the cushion must be properly sized to the individual.

### Limited warranty

Cushion: 24 months

Cover: 6 months



### Our evidence shows the ROHO HIGH PROFILE Single Compartment Cushion is ideal for people who:

- Currently have, have a history of or are at risk for skin/soft tissue breakdown and/or any category/stage pressure injury (including deep tissue injury)
- Require accommodation or correction of pelvic asymmetry, approx. 3 in. (7.5 cm)
- Need skin protection, maximum immersion and envelopment
- Have normal or impaired sensation
- Have a history or at risk of skin/soft tissue breakdown



## ROHO® MID PROFILE™ Single Compartment Cushion

The 3.25 in. (8.5 cm) individual air cells bring you the perfect balance immersion, protection and positioning control, redistributing pressure to decrease friction and shear. With a simple inflation adjustment, your ROHO conforms to your shape, minimizing tissue deformation and maximizing blood flow.



**Pressure injury prevention and treatment**



**Adjustable**



**Stability**



**Lightweight**



**Washable**

### Dimensions

Width: 13.25 in. to 20 in. / 33.5 cm to 50.5 cm

Depth: 13.25 in. to 20 in. / 33.5 cm to 50.5 cm

Height: 3.25 in. / 8.5 cm

Weight: Approx. 3.0 lbs. / 1.4 kg

### Weight limit

No weight limit.

However, the cushion must be properly sized to the individual.

### Limited warranty

Cushion: 24 months

Cover: 6 months



### Our evidence shows the MID PROFILE Single Compartment Cushion is ideal for people who:

- Currently have, have a history of, or are at risk for skin/soft tissue breakdown and/or any category/stage pressure injury (including deep tissue injury)
- Need skin protection, immersion and envelopment
- Have normal or impaired sensation



## ROHO® LOW PROFILE® Single Compartment Cushion

Designed for those with active lifestyles, the 2.25 in. (5.5 cm) tall, interconnected air cells maximize your mobility for transfers. The cushion provides the right balance of immersion and stability, protecting you from tissue breakdown by redistributing your weight across the cushion.

The LOW PROFILE cushion is also ideal for use as a secondary cushion on an office chair, in the car, or in the bath – anywhere you need extra cushioning and protection from injury.



**Lightweight**



**Transfer stability**



**Adjustable**

### Dimensions

Width: 11.5 in. to 23.5 in. /  
29 cm to 59.5 cm

Depth: 11.5 in. to 18.25 in. /  
29 cm to 46.5 cm

Height: 2.25 in. / 5.5 cm

Weight: Approx. 2.5 lbs. /  
1.1 kg

### Weight limit

No weight limit.  
However, the cushion  
must be properly sized to  
the individual.

### Limited warranty

Cushion: 24 months  
Cover: 6 months



**Our evidence shows the  
LOW PROFILE Single  
Compartment Cushion is  
ideal for people who:**

- Currently have, have a history of or are at risk for skin/soft tissue breakdown and/or any category/stage pressure injury (including deep tissue injury)
- Have normal or impaired sensation
- Have a history or are at risk of skin/soft tissue breakdown

# ROHO Cushion Sizing Chart

## Standard DRY FLOATATION® Cushion Sizing Chart for:

- ROHO LOW PROFILE  
Single or Dual Compartment Cushion
- ROHO MID PROFILE  
Single or Dual Compartment Cushion
- ROHO HIGH PROFILE  
Single or Dual Compartment Cushion
- ROHO ENHANCER Cushion
- ROHO MID PROFILE  
Single Compartment with  
Sensor Ready Technology
- ROHO HIGH PROFILE  
Single Compartment with  
Sensor Ready Technology

## ROHO SELECT Series Cushion Sizing Chart for:

- ROHO QUADTRO SELECT  
LOW PROFILE Cushion
- ROHO QUADTRO SELECT  
MID PROFILE Cushion
- ROHO QUADTRO SELECT  
HIGH PROFILE Cushion
- ROHO CONTOUR SELECT Cushion

## Select the correct size ROHO Cushion for your wheelchair:

### Step 1

Measure the wheelchair seat width and depth.

Note: Wheelchair manufacturers' published seat dimensions may vary.

### Step 2


Use your wheelchair seat dimensions and the below table to help you determine the number of cells in cushion width and in depth. The current version is available at ROHO.com.

## ROHO® Neoprene Cushion Sizing Chart

Fits wheelchair size (in. / cm)	Number of cells
10" / 22 to 26.5	5
11" to 12" / 26.5 to 31	6
13" / 31 to 35	7
14" to 15" / 35 to 39.5	8
16" to 17" / 39.5 to 44	9
18" to 19" / 44 to 48.5	10
20" to 21" / 48.5 to 52.5	11
22" / 52.5 to 57	12
23" to 24" / 57 to 61	13
25" / 61 to 65.5	14
26" to 27" / 65.5 to 70	15

Note: Availability of cushion sizes varies by model.

## Standard DRY FLOATATION Cushions

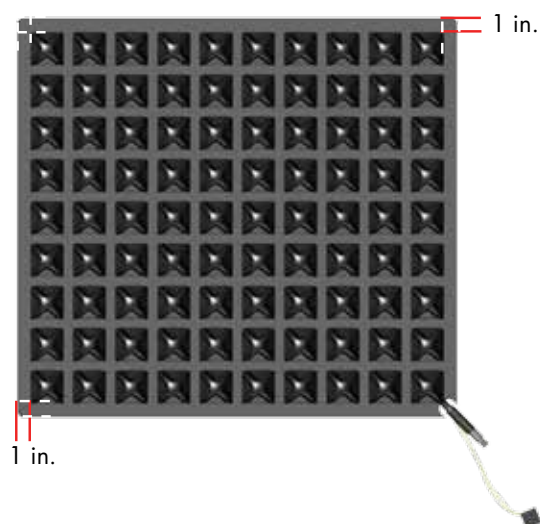

## SELECT Series Cushions

