



 Breezy

BREEZY ULTRA 4





Breezy® — the tradition of superior clinical benefit integrated with stylish looks continues with the new Breezy Ultra<sup>4</sup> high-strength, lightweight wheelchair. Compare these features and you will find the Breezy Ultra<sup>4</sup> sets a higher standard for performance, versatility, and value.

**Lighter Weight Aluminum Frame** — Every ounce of added weight requires greater energy and strength, not only to propel but also to load in and out of a car. Breezy Ultra<sup>4</sup>, made out of corrosion-resistant aluminum, is one of the lightest wheelchairs in its class.

**Quick-Release Wheels Are Standard** — An added cost on other comparable wheelchairs, all Breezy Ultra<sup>4</sup> users will appreciate the ability to remove the wheels with just the push of a button. This makes transferring into vehicles easy by converting the chair into three smaller, lightweight components.

**Enhanced User Function and Positioning** — The Breezy Ultra<sup>4</sup>'s package of adjustable features and options, plus the ability to accommodate Jay® seating products, allows you to create a custom fit to accommodate a variety of positioning needs. The Breezy Ultra<sup>4</sup> can be set up for foot propulsion, accommodate needs for added lateral stability, provide for planar flexion, and aid in positioning kyphotic postures making it one of the most versatile models available.



Colors



Actual colors may appear different



## Premium Features are Standard on Breezy Ultra<sup>4</sup>




### Quick-Release Wheels

Make loading in/out of vehicles easy.



### Angle-Adjustable Front Casters

Accommodate angled seating surfaces without moving the caster housings out of square.

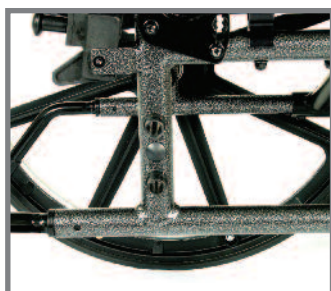

### Adjustable Axle Positions

Three seat height positions and the ability to add a fixed tilt allow for greater flexibility.



### Quick-Adjust Footrests

Adjust height without tools using spring-loaded detent buttons. Heel loops adjust with Velcro®.

## Key Options



### Angle-Adjustable Back

Back adjusts from 6° anterior to 24° posterior in 6° increments for better positioning. Height adjusts to 17", 18" or 19".



### Tension Adjustable Back

Highly adjustable velcro strap system provides customized postural support.



### Height-Adjustable Arms

4" vertical adjustment allows for use of cushions.



### Hemi Frame

Allows for a seat height as low as 13.5".

Specifications		
Frame Weight <sup>1</sup>	30.9 lbs (14.02 kg)	
Frame Material	T6 6061 Aluminum	
Seat Width & Depth	14"x16"	14"x18"
	16"x16"	16"x18"
	18"x16"	18"x18"
	20"x16"	20"x18"
Back Height	Adjustable	17", 18", 19"
Seat Height	13.5", 14.5", 15.5", 16.5", 17.5", 18.5", 19.5"	
Seat Options	Padded Seat Jay Basic Cushion	Non-Padded Seat
Back Options	Padded Back w/wo Pocket Tension Adjustable Back	Angle Adjustable Back Jay Basic® Back
Front Casters	6" Polyurethane 6" Solid	8" Polyurethane 8" Solid
Rear Wheels with Aluminum Handrims	20"/24" Full Profile Polyurethane 22"/24" Low Profile Polyurethane 24" Full Pneumatic w/ Airless Insert	
Axles	Quick Release	
Armrest	Flip-Back, Desk or Full Flip-Back Height Adjustable, Desk or Full	
Front Hangers	Swing-Away with Heel Loops Elevating Legrest	70° Swingaway Articulating Legrest
Maximum User Weight	250 lbs (113.4 kg)	
HCPCs Code	K0004 (meets requirements for)	
1 18"x16" Configuration without front hangers		
2 Only works on left, does not work on Angle Adjustable Back		

Accessories	
Seat Belt, Auto Style	BZU9001
IV Pole	BZU9003 <sup>2</sup>
Tray Table, Full Arm	BZU9002
Tray Table, Desk Arm	BZU9009
Oxygen Tank Holder	BZU9004 <sup>2</sup>
Crutch Holder	BZU9005

## Now providing BASIC seating comfort is Easy with Breezy.

The JAY® Basic and JAY® Basic Cushion can be ordered and shipped as part of a configured Breezy saving you time and money. Back hardware is pre-assembled onto the chair frame.



**JAY**  
Basic BACK

**Quick Ship or Built-to-Order** — The most popular configurations are stocked with front hangers in the box ready to ship allowing you to meet your customers' needs quickly with reduced inventory. We will also custom build a wheelchair with your choice of the available options.

**Custom Configurations:** Indicate your selection of standard features and options on the Breezy Ultra<sup>4</sup> order form.

**Quick-Ship:** These popular configurations are packaged in one box, ready for shipment from Sunrise Medical distribution centers.

## BREEZY ULTRA 4 SPECIFICATIONS

Quick Ship models include: 24" low-profile wheels with aluminum handrims, 8" polyurethane front casters, 17.5"-19.5" seat-to-floor height, nylon non-padded seat sling and back pocket.

### Breezy Ultra<sup>4</sup> 16"

B1616FAPDS	•	16"	16"	•	•	•	•	•	•	•
B1616HAPDSAT	•	16"	16"	•	•	•	•	•	•	•
B1616ZHAPDSAT	•	16"	16"	•	•	•	•	•	•	•
B1616ZHAPFSAT	•	16"	16"	•	•	•	•	•	•	•
B1616HAPFSAT	•	16"	16"	•	•	•	•	•	•	•

### Breezy Ultra<sup>4</sup> 18"

B1816FAPDS	•	18"	16"	•	•	•	•	•	•	•
B1816HAPDSAT	•	18"	16"	•	•	•	•	•	•	•
B1816HAPFSAT	•	18"	16"	•	•	•	•	•	•	•
B1816ZHAPDSAT	•	18"	16"	•	•	•	•	•	•	•
B1816ZHAPFSAT	•	18"	16"	•	•	•	•	•	•	•
B1818FAPDS	•	18"	18"	•	•	•	•	•	•	•
B1818HAPDSAT	•	18"	18"	•	•	•	•	•	•	•
B1818HAPFSAT	•	18"	18"	•	•	•	•	•	•	•
B1818HAPDSAT	•	18"	18"	•	•	•	•	•	•	•
B1818ZHAPFSAT	•	18"	18"	•	•	•	•	•	•	•

### Breezy Ultra<sup>4</sup> 20"

B2018ZHAPDSAT	•	20"	18"	•	•	•	•	•	•	•
B2018HAPDSAT	•	20"	18"	•	•	•	•	•	•	•
B2018ZHAPFSAT	•	20"	18"	•	•	•	•	•	•	•
B2018HAPFSAT	•	20"	18"	•	•	•	•	•	•	•

Seat-To-Floor-Height Matrix 6" Front Casters

Seat To Floor Height	Frame	Rear Wheel Size	Rear Wheel Axle Position	Front Caster Fork Position
13 1/2"	Hemi	20"	Top	Top- 6" Fork
14 1/2"	Hemi	20"	Bottom	Bottom- 6" Fork
15 1/2"	Standard	20"	Top	Top- 6" Fork
16 1/2"	Standard	20"	Middle	Bottom- 6" Fork
17 1/2"	Standard	20"	Bottom	Middle- 8" Fork
16 1/2"	Standard	22"	Top	Bottom- 6" Fork
17 1/2"	Standard	22"	Middle	Middle- 8" Fork
18 1/2"	Standard	22"	Bottom	Bottom- 8" Fork
17 1/2"	Standard	24"	Top	Middle- 8" Fork
18 1/2"	Standard	24"	Middle	Bottom- 8" Fork

Seat-To-Floor-Height Matrix 8" Front Casters

Seat To Floor Height	Frame	Rear Wheel Size	Rear Wheel Axle Position	Front Caster Fork Position
17 1/2"	Standard	20"	Bottom	Top- 8" Fork
17 1/2"	Standard	22"	Middle	Top- 8" Fork
18 1/2"	Standard	22"	Bottom	Middle- 8" Fork
17 1/2"	Standard	24"	Top	Top- 8" Fork
18 1/2"	Standard	24"	Middle	Middle- 8" Fork
19 1/2"	Standard	24"	Bottom	Bottom- 8" Fork



Customer Service: 800.333.4000  
www.SunriseMedical.com

©08.2011 | Sunrise Medical (US) LLC  
930854 Rev. H