



Self propelled

Easy to operate for the independent user.

Safe and simple transfer

The flat seat and the brake position facilitate lateral transfer.


Etac Clean 24"

This variant of Etac Clean is suitable for those who can operate the chair independently and want to maximise their mobility.

Easy to manoeuvre, for the independent user

The Etac Clean 24-inch wheels are positioned to make it easy to manoeuvre the chair without compromising user safety or comfort. It is easy to transfer to and from the chair thanks to the large, unrestricted seat in front of the wheels. The brake handles are also positioned below the seat so that they do not get in the way during transfer, while remaining easy to reach.

Etac Clean 24" has the same basic benefits as the other members of the Clean family – but allows the user to operate the chair independently.

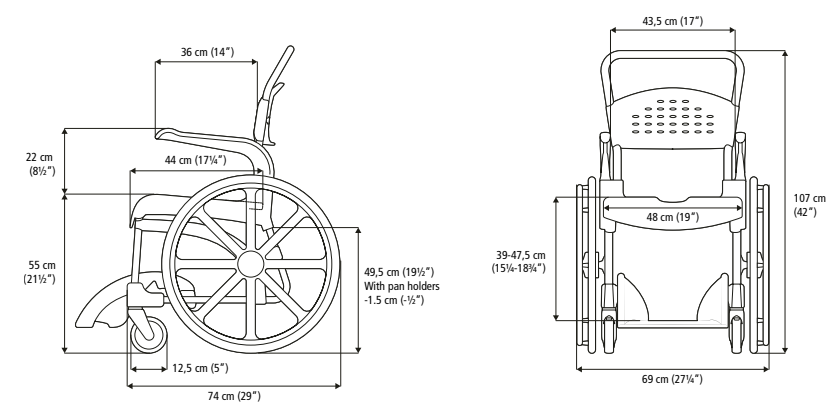

Etac Clean 24"

55 cm (21 1/2") with 24" rear wheels, white
55 cm (21 1/2") with 24" rear wheels, lagoon green

80229276
80229216

Weight
20 kg

Maximum user weight
130 kg

Material

Rear wheel adapter: aluminum
Rear wheels: polyamide, solid tires
Chassis: Powder-painted steel tubing
Seat/back support: Polypropylene
Arm supports: Polypropylene
Foot support: Polypropylene
Castor brackets: Polyamide
Castors, ball-bearing: Polyamide/TPE

Maintenance

Clean the product with a solvent-free cleaning agent (pH level 5–9), or with a 70% disinfectant solution. Can be decontaminated at max 85°C / 185°F if so required.



Accessories



Soft back support
80209225
Adjustable.
Washing at max 60° C (140°F).



Soft pad
Grey: 80209260
Heat-reflecting with a non-slip surface. Recess width: 17 cm (6¾"). Material: polyethylene. Cleaning max max. 60° C (140 °F) or with alcohol wipe.



Side Support
80209509
Adjustable in width. Material: Polyester fabric, Velcro, polyethylene. Washing at max. 60°C (140°F). *Can not be used in combination with: 80209225.*



Soft back pad
Grey: 80209261
Heat-reflecting with a non-slip surface. Material: polyethylene. Cleaning max. 60°C (140°F) or with alcohol wipe.



Crossbar
80209246
Attached to the arm supports. Material: Aluminium, cellular polyethylene, polypropylene. *Can not be used in combination with: 80209228, 80209237 and 80209226.*



Seat cover
80209247
When a complete seat is required. Material: Polyurethane. Cleaning max. 85° C (185°F) or with alcohol wipe.



Chest harness strap
80209442
Adjustable in length. Attach to soft back support. Material: Polyester.



Comfort Kit
80209237
Includes both Comfort Seat Soft and Comfort Cover. Cleaning max. 60° C (140 °F) or with alcohol wipe. *Can not be used in combination with: 80209225 and 80209246.*



Positioning belt
80209440
Adjustable in length. Attachment with velcro straps. Material: Polyester.



Comfort Cover
80209228
A soft back support with integrated arm support pads. Material: polyethylene. Cleaning max. 60° C (140 °F) or with alcohol wipe. *Can not be used in combination with: 80209225 and 80209246.*



Safety strap
80209034
Attachments on the rear side of the back support. Always use together with positioning belt. Material: Polyester. *Can not be used in combination with: 80209225.*



Comfort Arm Pads
80209226
Soft and warm, leaving the handles free to ensure a firm grip when standing up. Material: polyethylene. Cleaning max. 60° C (140 °F) or with alcohol wipe. *Can not be used in combination with: 80209246.*



Cushion, protect
80209245
Air-inflated ROHO cushion with excellent pressure distribution properties. Adjustable air pressure. Recess width: 21 cm (8¼"). Material: Neoprene. Cleaning max: 60° C (140°F).



Comfort Seat Soft
Small opening width 15 cm (6"): 80209301
narrow opening width 18 cm (7"): 80209227
A soft contoured seat. Material: polyurethane. Cleaning max. 85° C (185°F) or with alcohol wipe.

Accessories



Soap dish (does not fit Clean 24")
Grey: 80209266
Attaches to the side of the seat. Built-in shower handle holder. Material: Polypropylene. (Does not fit Clean 24")




Bucket holder
Grey: 80209258
Fits Bucket with handle and lid and other buckets with a thick rim. Material: Polypropene.



Pan
Grey: 80209254
Material: polypropylene.
Pan incl. lid with handle
Grey: 80209255
Material: Polypropylene.



Attachment for low mounting of pan/bucket
Grey: 80209259
Fitted to the pan/bucket holder. Lowers the pan/bucket holder with 6 cm (2½"). Material: Polypropene.

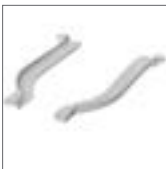

Pan holder
Grey: 80209257
Fits both pan, bedpan and regular hospital pans with a thin rim.



Disposable pans
80209282
Strong and firm. 50 pcs.



Bucket with handle and lid
80302026
The bucket holder is not included. Material: Polypropylene.



Disposable pan holder kit
Grey: 80209229
Fits Disposable pan.

Additional accessories – Clean Height Adjustable only

Additional accessories – Clean 49 and 55 cm only



Flip up foot supports
80209250
Removable, adjustable height. Attachments included. Material: powder coated aluminum and polypropylene foot supports.



Rear wheel adapter
White: 80209262
Material: Aluminium, polyamide.



Elevating leg support
Right: 80209252
Left: 80209253
Attachments included. Material: ABS plastic, stainless.



Heightening kit with wheels
White, 4 pcs: 80209511
Lagoon green, 4 pcs: 80209512
Raise the seat by 5 cm (2"). Material: Aluminium.



75 mm (3") castors
White, front, pair: 80209515
White, back pair: 80209516
Lagoon green, front, pair: 80209513
Lagoon green, back, pair: 80209514
Reduce seat height by 5 cm (2"). Mounted on Clean 44 cm as standard. Material: polyamide/TPE.