

SPECIFICATIONS

additional specifications at zingstanders.com

Weight Limit	Size 1- 70 lbs. (32kg) Size 2- 154 lbs. (70kg)
Height Range	Size 1- up to 44" (112cm) Size 2- 40"-60" (102-152cm)
Transfer Height	Size 1- 32"(81cm) Size 2- 31" (79cm)
Hip to Foot Plate Range	Size 1- 7"-22" (18-56cm) Size 2- 15"-33" (38-84cm)
Upper Body Support Height Range	Size 1- 4" (10cm) Size 2- 6" (15cm)
Knee Pad Depth Range	Size 1- 5" (13cm) Size 2- 7" (18cm)
Weight of Basic Unit	Size 1- 55 lbs (25kg) Size 2- 82 lbs (37kg)
Frame Footprint	Size 1- 20"x28" (51x71cm) Size 2- 26.75"x39.75" (68x101cm)
30° Independent Leg Abduction on each leg	Available on all models
Supine	0°-90° (not available on Vertical & Prone)
Prone	0°-20°-MPS, MPS TT; 0°-90°-Prone (not available on Supine, Supine TT, & Vertical)
HPCS	E0641-Multi-Position Stander MPS, MPS TT E0638-Single Position Stander Supine, Supine TT, Vertical, & Prone
Positions	Supine (size 1 & 2) Vertical (size 1) Prone (size 1 & 2)
Adjustable Tray	Available on all models
Swing-Away Tray	Available (not available on Prone and Vertical)
Support	Individual Positioning Components, Padded Strapping
Dual Control	Available (not available on Vertical)
Gas Spring Lockout	Available (not available on Vertical)
Transfer Aid	Flat-to-Load (not available on Vertical)
Tool Free Adjustment	Yes
Pow'r Up Lift Option	Size 2 Only


Frame Colors



GREEN



PURPLE

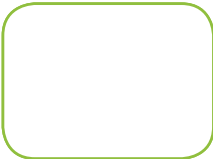

WHITE



ORANGE

Full Line of Standers

- MPS (multi-position stander)
- MPS TT (tilt table)
- Supine
- Supine TT
- Prone
- Vertical (size 1 only)



Zing MPS TT-Size 1

Hygienic Covers

- Available in 5 fun patterns for the pelvic, upper body, lateral, head and tray areas
- Prevents stander from getting soiled, especially during feeding
- Easily removed and washable



RAINBOW UNICORN



SEA CREATURES



FANCY PANTS ELEPHANTS



SPACE RACE



ALL SMILES ANIMALS

-MPS

-MPS TT

-Supine

-Supine TT

-Prone

-Vertical
(size 1 only)



Zing

zingstanders.com

Altimate Medical, Inc.
262 West First Street
Morton, MN 56270 USA
Tel: 800.342.8968
507.697.6393
Fax: 877.342.8968
507.697.6900
info@zingstanders.com
zingstanders.com




Visit our website at zingstanders.com for complete specifications on all Zing products and options.


Standing technology should only be used under the guidance of a physician with recommendations for standing program protocol and any medical precautions. Standing programs should be monitored by the attending therapist. AMI maintains a policy of continual product improvement and reserves the right to change features, specifications, and prices without prior notification. Check with AMI for latest information.
FORM ZINGMINICAT031419 REV A Copyright © 2019 Altimate Medical, Inc. All rights reserved. Printed in the USA.

Zing

Size 1
up to 44"
(112cm)

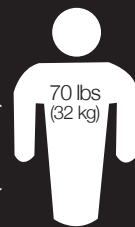

Size 2
40"-60"
(102-152cm)

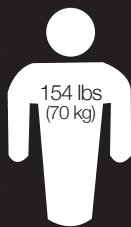

The Zing Advantage



Zing MPS-Size 2

True Multi-Position Stander

- The MPS and MPS TT allow supine, vertical and prone positioning without flipping foot plates or pads, moving the tray or transferring out and back in the stander
- Smooth transition from a flat-to-load transfer position, supine, vertical and 20 degrees prone - all in one motion thanks to a unique pivot
- Provides therapeutic opportunities not possible with other multi-position standers - use prone for head control and when fatigue sets in, simply move to a slight supine position to rest but remain in the stander
- Frame will only go supine and upright with Prone Position Lockout engaged

True Hip Pivot Leg Abduction

- Pivot points are in line with the hips, promoting proper alignment in the hips, knees and feet
- Feet maintain the angle with the rest of the legs - no wedges needed and no additional adjustment to the knees or foot plates
- No need for abduction pad/pommel due to the innovative design
- Up to 30 degrees of abduction per leg - accomplished through a single knob
- Promotes musculoskeletal development for children during critical growth years

- Abduction can be adjusted and used throughout the entire positioning range

Plenty of Growth

- Size 1: Fits children from infancy to 44", and up to 70lbs
- Size 2: Fits children from 40"-60", and up to 154lbs
- Can easily grow with the child and adapt to their changing needs
- Many supports, pads and straps are available in multiple sizes
- Comfortable transfer height of 32" makes transfers more manageable as children grow

Remarkably Modular

- Can be configured as a simple stander or fit a child with the most complex needs
- Designed to grow with the child - options can be added as needed
- Great for a multi-user facility or school



Zing Prone-Size 1



Zing MPS TT-Size 1

Multi-Use Tray

- The MPS, Supine, MPS TT and Supine TT feature a multi-use tray that is angle adjustable and swings away either direction
- Tray can remain attached during transfers
- More activities are possible by swinging tray away during standing therapy

Unique Lift Mechanisms

- Dual controls are available allowing angle adjustment by foot pedal or hand control
- Gas Spring Lockout prevents an inadvertent change of position, especially important when around other children
- The Pow'r Up Lift option on Size 2 brings the user to the standing position with the touch of a button

Accommodating Contractures

- No special modifications are needed to accommodate children with contractures
- Using the multi-adjustable foot plates, multi-adjustable knee pads, and the independent calf pads, it is possible to adjust for contractures in one or both sides

Easy Adjustment

- Tool-free adjustment of all components
- Every adjustment can be made with a simple turn of a knob or ratchet handle



Zing MPS TT-Size 2

SUPINE TO PRONE

MPS Models Size 1 and Size 2



Supine Lying - Flat-to-Load



Supine



Vertical



Prone



Zing MPS-Size 1