

SPECIFICATIONS

<b>Model Number</b>	<b>500</b> stretcher/crib
Patient Surface	26" x 56¼" (66 x 143 cm)
Overall Length	60" (152 cm)
Overall Width	30¼" (77 cm)
<b>Height</b>	
Fixed Litter Top	33¼" (84 cm)
Top of Crib	59½" (151 cm)
With Security Top	76¾" (195 cm)
<b>Litter Positioning</b>	
Fowler Angles	0°, 20°, 30°
Trendelenburg	0°, 5°, 10°
Fowler Angles	-5°
Side Rails Intermediate Stop	13¾" (35 cm)
Casters	5" (13 cm)
Warranty	1 year parts and labor

- 200 lbs. weight capacity rating
- Two 5" locking casters and one steering caster for increased maneuverability
- Trendelenburg/Reverse Trendelenburg
- One handed rail operation
- White powder coat finish
- Unique rounded corners
- Rails with 2" spacing between bars
- Two IV pole receptacles
- Full function fowler
- Litter Top: Solid pan
- Color decals: Blue

✓ **LATEX-FREE**



OPTIONS

Accessory	Accessory
<b>Utility Shelf:</b>	<b>IV Pole:</b>
Utility shelf .....5103001	IV Pole .....5116001
<b>Security Top</b>	<b>Mattress and Pad:</b>
Windowed top with 6 mounting clips .....5107001 <i>(will accommodate a mist tent - requires 5104002 for mounting)</i>	Crib side rail bumper pad (pair) .....5109002
Clear top with 6 mounting clips.....5117001 <i>(will accommodate a mist tent)</i>	Pressure Care Mattress 4" .....5817002
<b>Adapters</b>	• <i>Mattress meets fire codes TB-117.</i>
Traction and croupette adapter.....5104002	Mattress with fire barrier layer 3½" ..... 5860003
Traction adapter, pair .....5119001 <i>(includes mounting kit)</i>	• <i>Mattress meets fire codes TB-117 &amp; TB-129.</i>
Croupette humidifier adapter .....5108001	• <i>For additional fire layer barrier protection, contact customer service for a complete listing of fire codes met by this mattress.</i>
Infusion pump adapter.....5106001	Mattress 3½" (standard) .....5860002
<b>O<sub>2</sub> Holder:</b>	• <i>Phthalate free • compliant with HR 4040.</i>
Oxygen tank holder .....5114001	<b>Color Decals:</b>
Oxygen tank holder (for larger regulators).....5114001-1	Red Head/Foot End Labels .....5118001
	Blue Head/Foot End Labels (standard) ....5118002
	Green Head/Foot End Labels .....5118003



# stretcher series 500

## Function as Appealing as its Form

Pedigo's 500 Pediatric Crib Stretcher, with its full function capabilities represents the standard of excellence in hospital cribs. The aesthetically pleasing design makes it perfect for any facility.



### Casters

As a pediatric transport stretcher, the 500 rolls easily on 5" oversized casters. Two of the four casters lock to secure the stretcher while a third steering caster may be engaged for increased control and maneuverability.



### Total Accessibility

The unique rounded corners reduce the overall size of the 500, yet provide the same amount of useful area on the 26½" x 55¾" mattress.

The two-post support system provides 97% access to the patient by the medical team.



### Side Rails

Side rails have two-inch bar spacing and a special epoxy finish that is warmer to the touch than chrome.



### A True One-Handed Crib

For complete containment of the child, simply raise the counter balanced side rails to enclose the crib. Because the rails are counter balanced, there will be no sudden drop when released, preventing injury to the child or provider.

Dual-action release handle is childproof and easy to operate, even with a child in one arm.



Crib Security Handle



Squeeze Handle Together



Lift Rail Upward or Downward

**LATEX-FREE**

## stretcher options




### FIXED SECURITY TOP

Clear molded security top has an exclusive quick disconnect for fast, easy removal.



### SECURITY TOP - REMOVABLE

Clear molded security top has removable window for access into the croupette and crib.

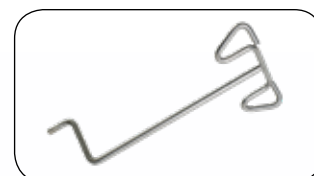

### IV POLE

IV poles allow for bottles or bags to be suspended inside or outside the stretcher area.



### UTILITY SHELF

Located at the base of the stretcher, utility shelf may be used for storing blankets or the child's personal items.



### O<sub>2</sub> HOLDERS

O<sub>2</sub> holder for E-size oxygen tank.

Also available is Pedigo's 5114001-1 O<sub>2</sub> Tank Holder to accommodate tanks with longer O<sub>2</sub> regulators.



### CRIB BUMPER PADS

Siderail Bumper Pads measures 2" thick and 7" high. Standard, baby soft white color.

