



## Viper Plus GT Wheelchair with Universal Armrests

HCPCS: K0004-1, K0004+K0195+E0993



FIGURE A

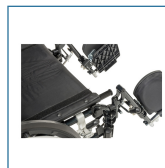


FIGURE B



FIGURE C



FIGURE D



FIGURE E



FIGURE F



FIGURE G



FIGURE H

- The new Viper Plus GT with Universal Armrests allows the user to easily convert the armrests into desk arm, half arm, and full arm with the push of a button. The armrests also adjust to eight individual heights in  $\frac{1}{2}$  increments for complete customization.
- Adjustable height and length, flip-back, removable, padded arms
- Built-in seat rail extensions and extendable upholstery easily adjusts seat depth from 16" to 18"
- All-Aluminum frame
- Dual Axle allows for multiple seat-to-floor height positions
- Over-center cross brace eliminates need for seat guides Works with seating and positioning accessories
- Adjustable angle back easily adjusts from 5° - 20°
- Adjustable height back adjusts from 17" - 19" in 1" increments, adaptable for individual patient needs
- Durable, heavy-gauge flame-retardant nylon upholstery with hook and loop fastener flaps
- Quick release, 24" composite Mag-style wheels are lightweight and maintenance free
- Adjustable angle caster forks have three height adjustments and angle adjustability
- Tool-free adjustable swing-away footrests also feature easy-to-access release handles
- Urethane rear tires and casters offer superior performance, smoother ride and are lightweight
- Adjustable anti-tipper with flip-up wheels Never needs to be removed
- Comes with push-to-lock wheel locks

### PRODUCTS

### Description

PLA416FBUARAD-ELR	ViperPlsGT,16"FlpBk.UniArm,ELR,1/ea
PLA416FBUARAD-SF	ViperPlsGT,16"FlpBk.UniArm,SF,1/ea
PLA418FBUARAD-ELR	ViperPlsGT,18"FlpBk.UniArm,ELR,1/ea
PLA418FBUARAD-SF	ViperPlsGT,18"FlpBk.UniArm,SF,1/ea
PLA420FBUARAD-ELR	ViperPlsGT,20"FlpBk.UniArm,ELR,1/ea
PLA420FBUARAD-SF	ViperPlsGT,20"FlpBk.UniArm,SF,1/ea
PLA422FBUARAD-ELR	ViperPlsGT,22"FlpBk.UniArm,ELR,1/ea
PLA422FBUARAD-SF	ViperPlsGT,22"FlpBk.UniArm,SF,1/ea

### SPECIFICATIONS

### Description

Arm Height From Floor	28.5"-31.5"
Anti-Tip Wheels	2"
Arm Height From Seat	8.5"-11.5"
Armrest Length	13.3"
Back Angle	5°-20°
Back of Chair Height	17" - 19"
Brakes	Push to lock

SPECIFICATIONS	Description
Color	Deep Green
Folded Dimensions	33"x12.25"x36"
Front Caster Diameter	8"
Handle Height Range	36"-38"
Material	Aluminum
Overall Length with Riggings	46"
Product Weight w/out Footrests	Seat Width 16": 34.5 lbs   Seat Width 18": 35.25 lbs   Seat Width 20" and 22": 36 lbs
Rear Wheels	24"
Seat Depth	16"-18"
Seat To Floor Height	17.5" - 19.5"
Seat Width	16"   18"   20"   22"
Warranty	Limited Lifetime
Weight Capacity	16", 18" and 20": 300 lbs   22": 350 lbs

ACCESSORIES	Description
STDS850	Seat Belt, Clasp Type,1/ea, EA
STDS851	Seat Belt, Velcro Close,1/ea, EA
STDS801GT	WheelLockExtension,ViperGT,SSP2,1/ea, EA
WASR	Limb Support, Right,1/cs, CS
WASL	Limb Support,Left,1/cs, CS
STDS840	RearWheel,TreadedUrethane, Viper GT 1/ea, EA
STDS1034	Wheelchair Cane/Crutch Holder,1/ea, EA
STDS828GT	AmputeeAdaptersViperGT,1pr, PR
STDS4056	One Arm drive for ViperGT, EA
STDS860	Forwrd Stablzrs/Univ.Fit ,1pr/bx, BX
STDS820	Telescoping IV Pole,W/C, 1/bx, BX
STDS834	OverheadAnti-Theft Bar,Single,1/ea, EA
STDSSHK	Super Hemi Kit Viper / Viper GT, 1kit/cs, CS
STDS807	Anti-Tips,Adj.w/Whls,Blk,1pr/bx, BX