



Power wheelchairs
M-Series




permobil

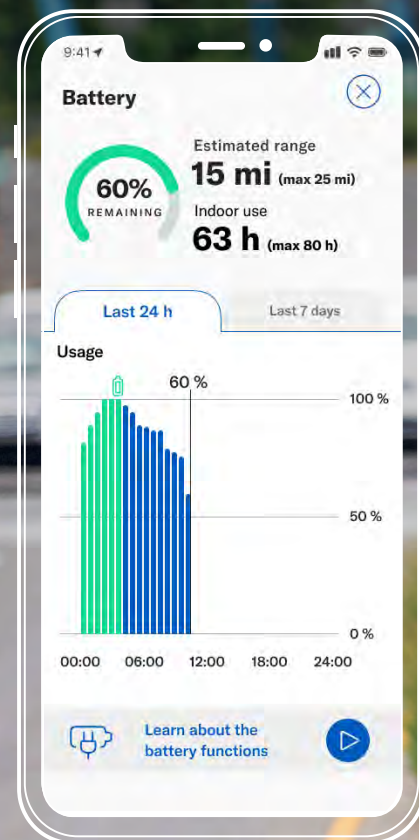

Battery status and travel range estimates help you know how far you can go

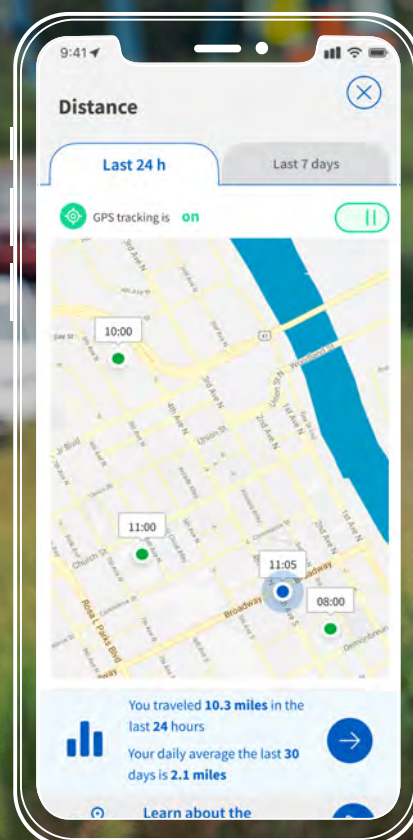


Integrated map illustrates daily and historic travel behavior

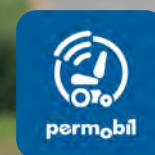

MyPermobil

MyPermobil connectivity is an available option on all new M-Series models.

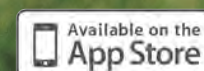

More adventure

The new Permobil battery and VoltPro charger work together to provide you with up to 25% greater driving distance, reduced charging times, and longer battery life. Use the MyPermobil app to get even more from every adventure with seat activity tracking, battery range estimates, and an integrated map.



Driving range up to 20 miles



Reduced charging times



Longer battery life

More freedom

Corpus®, our most capable power wheelchair seating system, now provides you with even more positioning options. Enhanced stability on the new M-Series enables you to optimize your repositioning, functional reach, and power drive positions for even greater freedom of movement.

Up to 3.2 mph
elevated drive
speed

Full suspension at any
elevation, seat position,
or drive speed





180° of power recline
50° of posterior tilt

90° – 180°
of legrest elevation

9:41

Seating

To save a seated or relaxed position, move your seat to the desired position and press the corresponding button. The slider will then indicate which position you are currently closer to.

Your current position

Position 1

Position 2

Tilt

Pos 1: 5°

Pos 2: 35°

35°

Recline

Pos 1: 95°

Pos 2: 115°

115°

Legrest

Pos 1: 95°

Pos 2: 135°

135°

Elevation

Pos 1: 1"

Pos 2: 0"

0"

88 seat adjustments today

Last 7 days average is 68

→

Learn about the seat functions

▶



MyPermobil

Live seating angles
and activity tracking

Reach for more

ActiveHeight™ and ActiveReach™ combine seamlessly to give you the ultimate advantage in accessing your environment. Seat elevation with integrated forward tilt creates the most innovative, functional, and stable seating solution on the market today.

	M3	M5
↑ ActiveHeight	12"	14"
↶ ActiveReach	20°	20°

10° of no charge ActiveReach with ActiveHeight



3.5° rearward shift during elevate for added stability and reduced footprint

12"

3 points of contact for added stability

Conquer more

FlexLink™ and FlexLink Pro™ suspension balances stability, traction, and ride comfort for optimal performance and advanced vibration reduction. Factory optimized and adjustable suspension adapts to your needs and environment. Be ready to handle any adventure with high-powered LED lights to guide your way.



FlexLink
M3 Corpus



FlexLink Pro
M5 Corpus

2X more powerful
LED lights for true
pathway lighting and
enhanced visibility



M3 Corpus

Meticulously engineered to help you maintain better positioning, feel more secure, and reduce driving fatigue so you can conquer daily activities with ease. Its compact base offers maneuverability in tight spaces without sacrificing performance. Featuring the innovative FlexLink suspension geometry with independent caster swing-arms helping to keep you stable and in control as you traverse the obstacles in your world.



M5 Corpus

Take control with high-performance motors enabling a top speed of 7.5 mph and exceptional low-end torque that lets you pivot like a pro. Experience unsurpassed stability and ride comfort with FlexLink Pro independent drive wheel suspension. The M5 Corpus will take you to new heights with increased stability and full suspension even in 14" of seat elevation.



BodiLink on Corpus


Specially designed for the Corpus seating system, BodiLink accessories feature best-in-class hardware. Offering true anatomical support with the highest level of adjustability and comfort to meet your individual positioning needs.

Head Support—
Easy to adjust piano hinge and double-taper joints for maximum strength


Lateral Trunk Support—
Designed for enhanced comfort and trunk alignment with unique wrist joint rotation adjustability allowing pad to match body contours

Lateral Pelvic/Thigh Support—
Easily position along the UniTrack rail with 360° pad adjustment and user-friendly quick release lever

Accent options

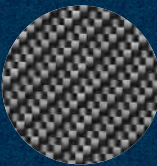

Carbon Fiber

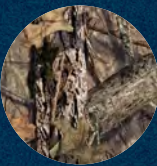

Mossy Oak® Breakup® Country™



Water Lily



Fireball



CDG Poly

Midnight

Arctic

Violet

Sky

Aqua

Spring

Canary

Sunset

Chili

Flamingo

Platinum

Bronze

Gold

Please note: Actual product colors may vary from colors shown



Redesigned essentials carrier provides the ultimate storage solution

Compact shell and smooth interior maximizes storage space while minimizing overall length

Specifications

	M3 Corpus	M5 Corpus
Maximum user weight	330 lb	330 lb
Maximum speed	6 mph	7.5 mph
Range	16 mi (Grp 34) 20 mi (Grp 24)*	20 mi*
Base width	24"	25"
Base length	35.5"	37.5"
Base turning radius	20.5"	21.5"
Ground clearance	3"	3"
Recommended battery type	Group 34 or 24	Group 24
Weight including batteries	385 lb (Grp 34) 405 lb (Grp 24)	423 lb
Drive electronics	R-net 120A	R-net 120A
Seat-to-floor height	17.5", 18.5", 19.5"	17.5", 18.5", 19.5"
Seat elevation	Optional 12"	Optional 14"
Drive wheel	14"	14"
Casters	7"	8"
Tilt options	Posterior: 0° – 50° ActiveReach: 5°, 10° & 20° (up to 300 lb)	Posterior: 0° – 50° ActiveReach: 5°, 10° & 20° (up to 300 lb)
Legrest angle - power	85° – 180°	85° – 180°
Backrest angle - power recline	85° – 180°	85° – 180°
Backrest height	20", 23" – 28" (1" incr.)	20", 23" – 28" (1" incr.)
Seat depth	14" – 22" (1" incr.)	14" – 22" (1" incr.)
Seat width	17" – 23" (2" incr.)	17" – 23" (2" incr.)

* Actual driving range will vary based on driving and battery conditions.

Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and include base, batteries, power seat functions, and seating system without accessories.