


Adapt 2


Modular system

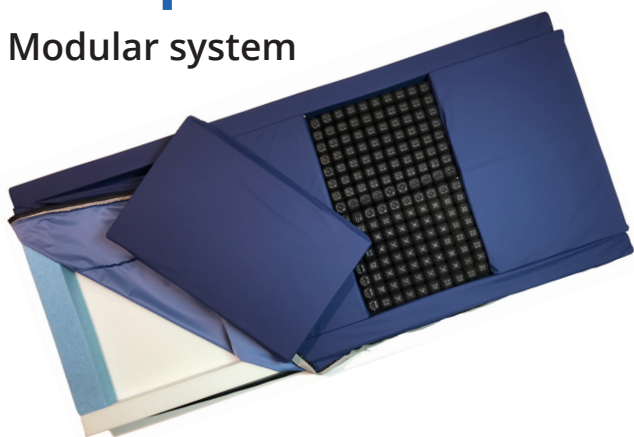


At a glance

- Ideal for long-term care clients
- Modular system uses ROHO mattress sections and foam inserts in a variety of configurations
- Wide transfer edges to help reduce the risk of falls
- Shear Reduction Technology (SRT)
- Full enclosure cover with Recovery 5 top
- Maximum client weight of up to 350-500 lb
- 3 widths: 36" / 39" / 42" x 80"L x 6"H
- 24-month warranty


The Adapt 2 mattress with SRT is an innovative way to reduce the effects of pressure, friction, and shear within the body of the mattress. The 3" wide firm perimeter borders and 4" permanent gel foam foot border hold the ROHO mattress sections or foam inserts in place without the need for snaps. The surface inserts can be configured based on the individual client's needs - see below for sample system configurations.

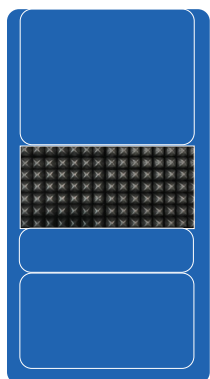


System 1

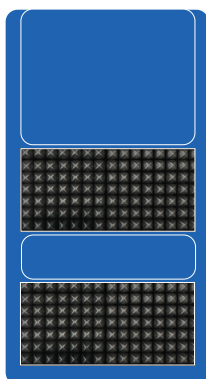


System 2

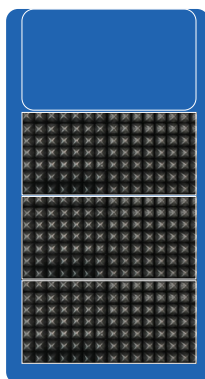

System 3

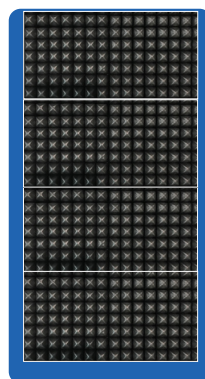

System 4

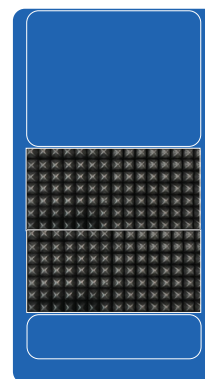

System 5