




Permobil, Atlas Hybrid Therapeutic Surface, ATL3580

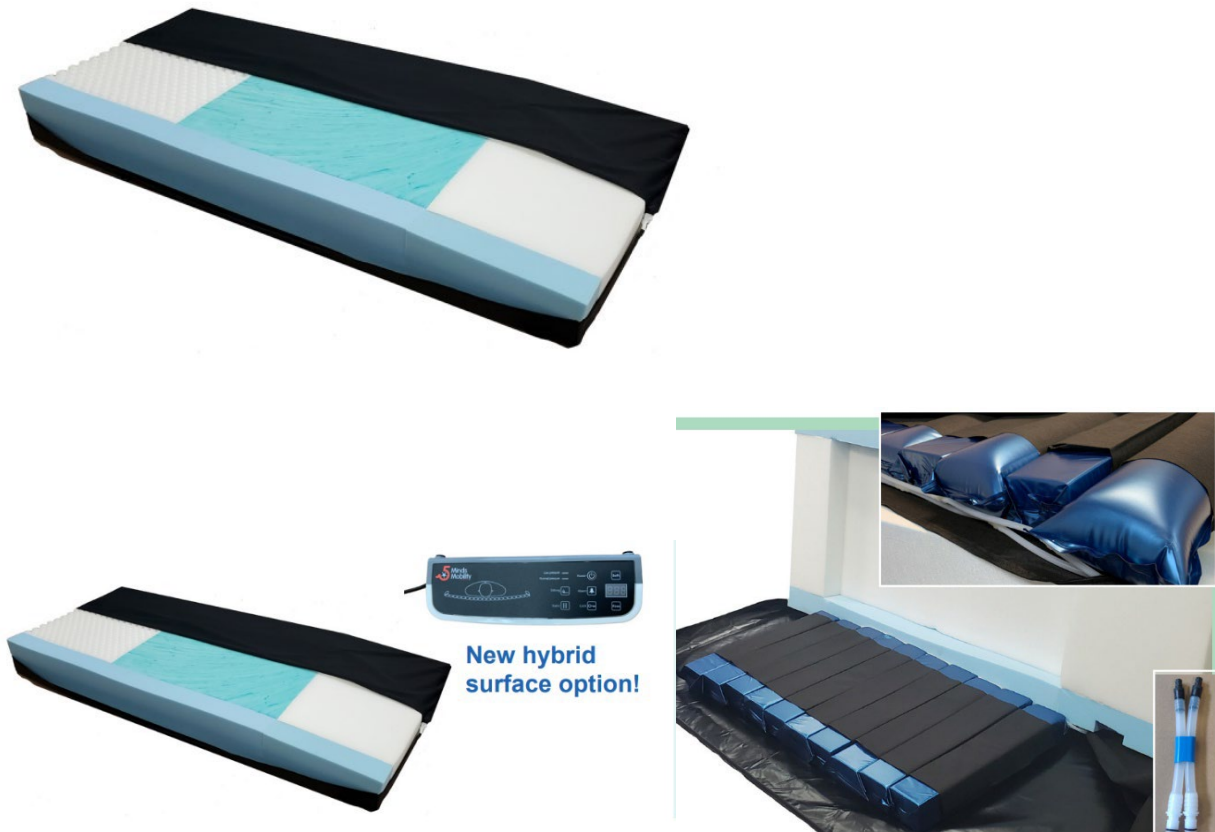

The Atlas is a non-powered, self-adjusting mattress that utilizes equalization technology to provide maximum pressure redistribution for the prevention and treatment of pressure injuries.

Starting at the head, the convoluted high resiliency foam improves client upper body mobility. Next is a cool climate visco gel swirl foam top layer at the torso and middle section for immersion and heat dissipation, followed by a 5-degree sloped foot section made with Serene™ foam to help transfer excess pressure from the heels. Beneath the top layer, interconnected chambers provide pressure redistribution, and there are 4" wide firm borders to ensure a solid transfer edge.

Atlas Hybrid Therapeutic Surface features:

- Suitable for vulnerable clients at risk for pressure injuries
- Interconnected chambers and valves for pressure redistribution
- Converts to a powered air surface with APT
- Sloped foot section helps reduce pressure from heels
- 4" wide firm borders for safe transfers
- Startex fluid-resistant, four-way stretch cover

Specifications:

- Maximum weight capacity : 700 lb
- Atlas is 7"H, available in six widths : 35" / 39" / 42" / 48" / 54" / 60"
- Two lengths : 80" / 84"