Plastic Transfer Bench

EN

REF	Description	Max User Weight
081295252	Plastic Transfer Bench	400 lbs / 181 kg

Before Assembly

Remove contents from package, making sure that all components are included and that no parts have been damaged during handling and shipping..

Familiarize yourself with the parts. Read the following instructions and make sure you understand them. For clarification of any step, consult a caregiver or healthcare professional.

Assembly Instructions

1. Place top of transfer bench on an even surface. Determine the side of bench that will go into the bathtub, and on that side attach the legs that have large suction cups on the bottom. Use the push pins, depressing them and slipping leg into the frame connector on the bottom side of the bench. Follow the same steps to attach the other suction-cupped leg (figure 1).
2. On the opposite side of bench, attach the legs that have small rubber tips. Depress push pins and insert one leg at a time into each of the corresponding receptacles (figure 2). Make sure the legs are secure before turning over the bench.
3. The armrest should be attached on the same side of large suction-cup legs. Hold the armrest so that it curves away from the bench, and then connect it to bench by depressing the push pins, inserting the arm and releasing the pins (figure 3).
4. Before attaching the backrest to the bench, assemble it using the push-pin method. Push the pins and slip the poles into the backrest poles (figure 4). Depress the push pins again to slide the poles of backrest into the holders underneath the seat (figure 5). To ensure security the backrest must go through both holders. Reverse the back to change entry from right to left or vice versa.
5. Position large suction-cup legs in the bathtub and adjust for individual user. To remove or reposition the legs, pull the tab on the suction cups. Check that the bench is level to be on the safe side.

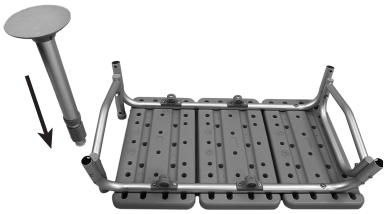


Figure 1

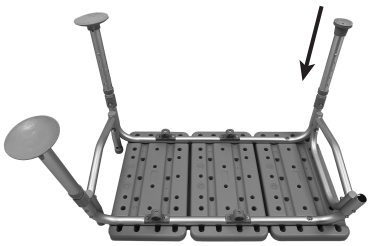


Figure 2

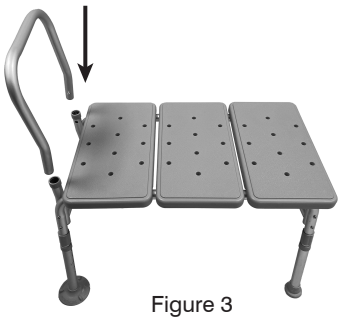


Figure 3

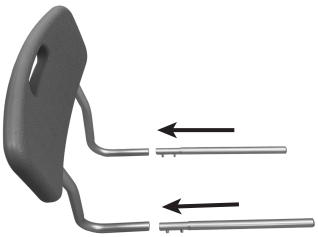


Figure 4

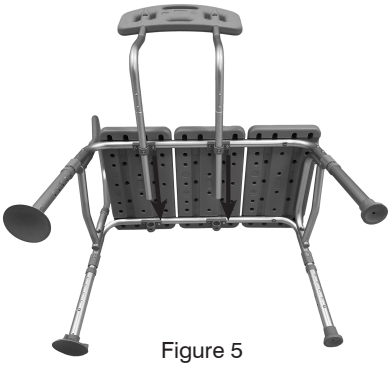


Figure 5



Figure 6

Care

Use a mild cleaning product—no abrasive, scouring pads or solvents including lacquer thinner, turpentine or acetone.

⚠ Precautions

Do not exceed the user weight capacity of 400 lbs (181kg). Test for stability before use.

Ensure that all push buttons are locked in place before use.

Limited Product Warranty

Performance Health Supply, LLC warrants that all products shall be free from manufacturing defects for a period of one year from the date of shipment. OTHER THAN THE FOREGOING WARRANTY, PERFORMANCE HEALTH SUPPLY, LLC MAKES NO REPRESENTATION OR WARRANTY OF ANY KIND, EXPRESS OR IMPLIED, AS TO MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE OR ANY OTHER MATTER RELATING TO THE PRODUCT. Customer must notify Performance Health Supply, LLC of any manufacturing defect within five days of discovery of the defect. Upon a determination by Performance Health Supply, LLC that the product is defective, Performance Health Supply, LLC sole and exclusive obligation shall be to provide customer with a replacement product or a full refund of the purchase price, at the sole discretion of Performance Health Supply, LLC, and to reimburse customer for any return freight costs. In no event shall Performance Health Supply, LLC be liable for any other costs or expenses including lost profits or direct, indirect, consequential, incidental or punitive damages. The warranty set forth herein shall govern all purchases, notwithstanding any inconsistent terms and conditions of customer's purchase order.

Serious events experienced from using this product should be reported immediately to Performance Health and to the local Competent Authority.

⚠ WARNING: This product can expose you to chemicals including Di(2-ethylhexyl) phthalate (DEHP) which is known to the State of California to cause cancer and or other reproductive harm. For more information go to www.P65Warnings.ca.gov

Distributed in the USA By:
Performance Health Supply LLC
W68 N158 Evergreen Blvd
Cedarburg, WI 53012, USA
+1800-323-5547
www.performancehealth.com

Distributed in Canada By:
Performance Health Canada Inc
6715 Millcreek Drive, Unit 2
Mississauga, ON L5N 5V2 Canada
+1800-665-9200
www.performancehealth.ca

PERFORMANCE HEALTH®