

ARJOHUNTLEIGH

GETINGE GROUP

VOYAGER PORTABLE CEILING LIFT



...with people in mind

PORTABLE CEILING LIFT VOYAGER


The *Voyager*® portable ceiling lift combines a solidly engineered design with features you would never expect to find in a portable lift.




Our unique Quick Release feature allows the user to disengage the strap and quickly install or remove the lift from the track.



The handset offers a degree of protection against splashes (IPX4)

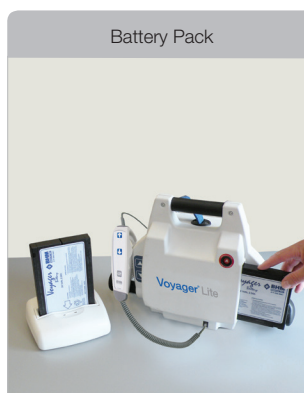

A 24-volt battery allows you to perform all your tasks worry-free. When necessary, it recharges in less than 2 hours.

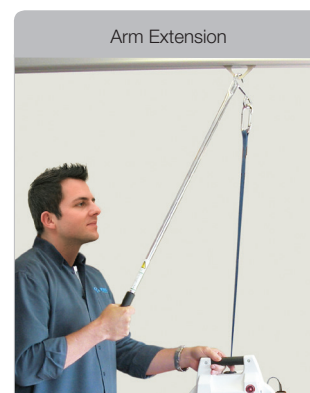

Combined with the Quick Release function, the arm extension gives the user the ability to detach the lift from the track. At the same time, unnecessary strain is avoided for the caregiver.

Compatible Slings – ArjoHuntleigh Ergofit slings are designed to be easily installed and removed.



Hammock Sling
THA model



Quick-Fit Sling
TIR model



Hygienic Sling
THY model



Walking Sling
TEM model

Product Information	
Weight	12 lb (5.5 kg)
Lifting Capacity	440 lb (200 kg)
Soft-start & soft-stop motor control	
Emergency Lowering Device	
Emergency Stop	
Current Limiter	
Vertical Display Speed	1.7 in/sec (4.3 cm/sec)
Up to 30 Transfers with Load of 180 lb (80 kg) on a Fully Charged Battery	
Strap Length	88" (2.2 m)
Automatically Stops if the Strap is Twisted	
Aluminium Frame	
Fire Retardant Plastic Cover	

Use only ArjoHuntleigh designed parts, i.e. parts specifically designed for the purpose, on equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development, we reserve the right to modify designs and specifications without prior notice.

© and ™ are trademarks belonging to the ArjoHuntleigh group of companies. © ArjoHuntleigh, 2014