

M1

Quality meets value.



permobil

Permobil tough

Finally, armrest and legrests designed to withstand the rigors of daily repositioning and transfers. Unlike traditional, wobbly rehab arms, the M1 locking flip back armrests are durable and secure. The center mount and multi-angle adjustable swing-away legrests are built to last with a full range of options to address your lower extremity positioning needs.

1" inch tubular frame
compatible with
virtually any seating



Locking flip back arms with
easy to release mechanism




ROHO® and Comfort
available on-chair



Swing in/out legs with
two-function release



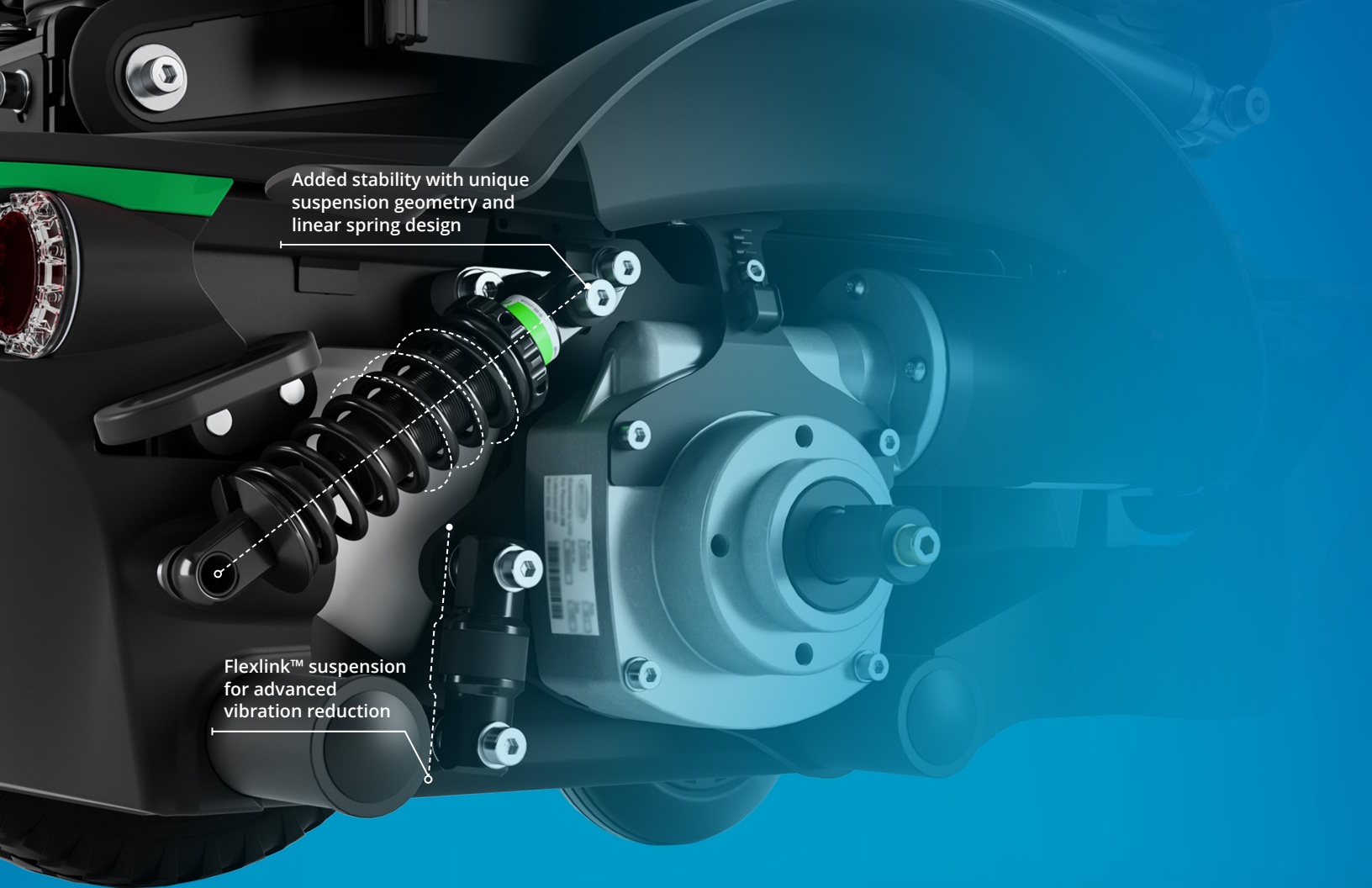


Added stability with unique suspension geometry and linear spring design

Flexlink™ suspension for advanced vibration reduction

Proven technology

Built on the modular chassis of the Permobil M3 Corpus®, the M1 features the same reliable drive units, patented FlexLink™ suspension and precision engineered CasterGlide™ technology to help ensure smooth, confident operation in a variety of driving scenarios.

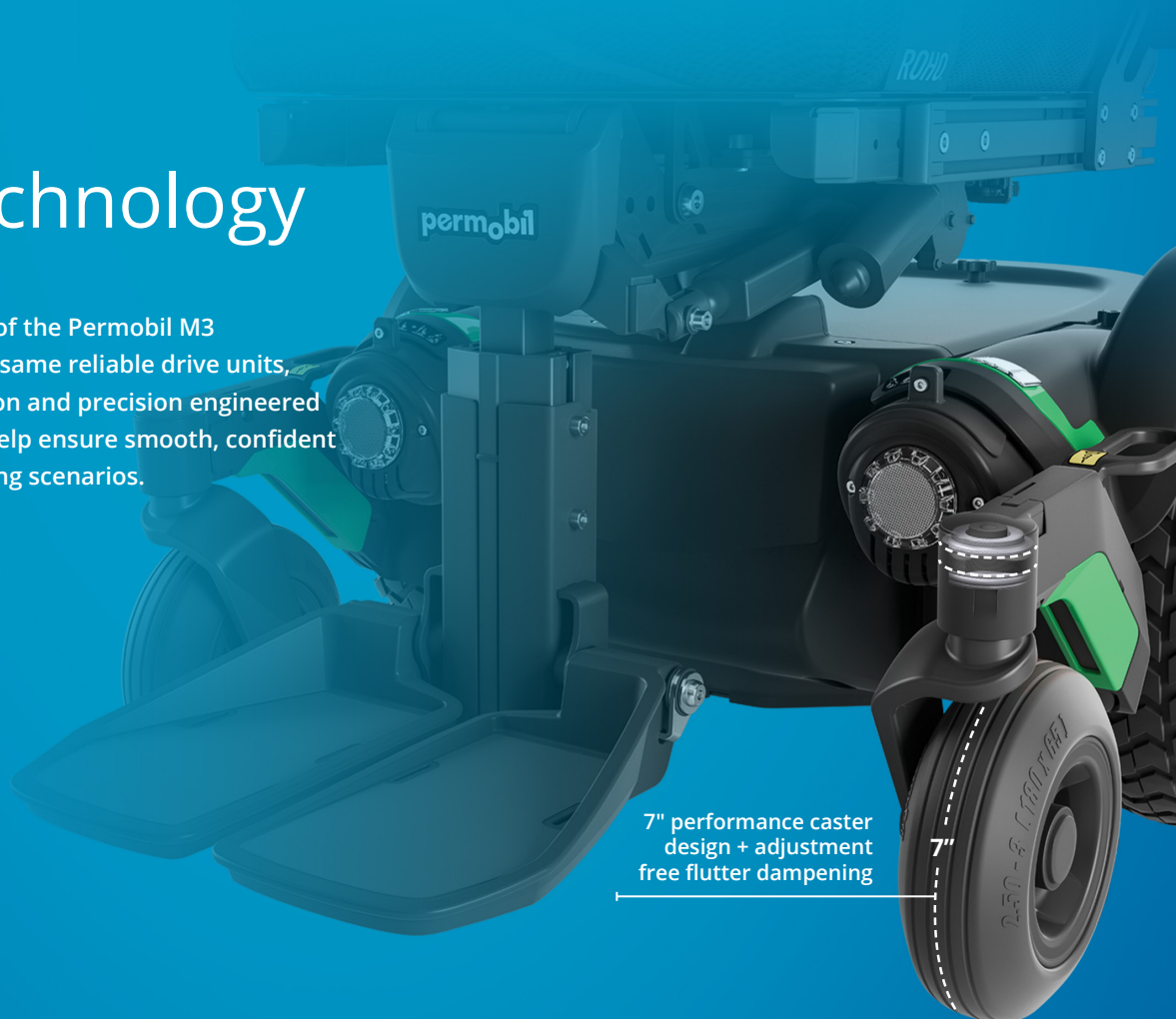

7" performance caster design + adjustment free flutter dampening

7"

Ready to go

Transport friendly features including standard crash tested transport brackets*, optional EZ fold backrest and a lower overall chair weight of 319 lbs with batteries, expands available travel options.

**Crash Tested to ISO 7176-19:2008(E)*



Specifications

Maximum user weight	300 lbs ^{††}
Maximum speed	5 mph
Range	16 miles*
Base width	24"
Base length	35.5"
Minimum turning radius	20.5"
Ground clearance	3"
Recommended battery type	Group 34
Weight including batteries	319 lbs**
Drive electronics	R-Net 90A, VR2 PM 90 (Advanced) or PM 70 (Advanced and Basic)
Seat-to-floor height	16", 17", 18"
Tilt options	0° – 50°

Recline options	85° – 110° (manual)
Legrest angle	Center Mount 80° – 140° (power) [>19" seat depth] 90° – 150° (power) [≤19" seat depth] 80° – 135° (manual) Swing-away 65° – 80°
Backrest height	18" – 29" (by 1" incr.)
Seat depth	14" – 22" (by 1" incr.)
Seat width	14" – 22" (by 2" incr.)
Armpad options	2"x10", 2"x14", 4"x10", 4"x13"

[†]Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions.

^{††}The maximum user weight for a seat depth of 14" cannot exceed 165 lbs.

**Actual driving range will vary based on driving and battery conditions.*

***Includes weight of base, curved backrest and batteries.*

Color options



Radiant Red



Sunburst Orange



Power Pink



Galactic Green



Ocean Blue



Midnight Black



Carbon Fiber



CDG Poly



Mossy Oak® Breakup® Country™



American Liberty



Midnight Patriot

Hydrographic options