WESTECH HEALTH CARE LTD.



Redefining Secure, Safe and Clean



SNAP 2 SECURE



Clean Shield elevated toilet seat



Support Arms for Clean Shield



Hinges and Rings



4.5" Raised Toilet Seat with Dual Lock Security™




Personal Wash Bidet™

SNAP 2 SECURE™

Your Security is Our Priority!

Raised toilet seats must be securely installed on the toilet to be safe and effective for the user. Our patented Snap 2 Secure installation system is engineered to secure the seat so it will not loosen with use, ever. It is so effective, uneven pressure can be applied to the arms and the seat will not loosen. All Bemis Independence seats have the Snap 2 Secure system.

How it works

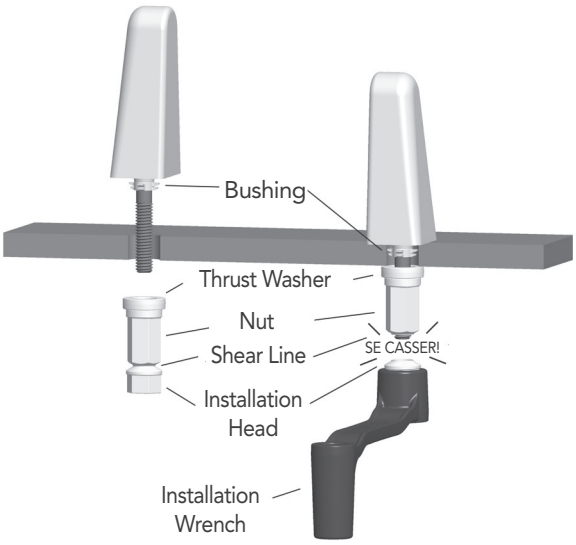

Bushing: Prevents the bolt from wiggling
Thrust Washer: Maintains the clamping force
Nut: Holds the seat in place
Shear Line: Installation Head snaps off here
Installation Head: Snaps off when the appropriate clamping force is reached
Installation Tool: Supplied with Clean Shield or use a 5/8" wrench



Hand Tighten



Snap off the Installation Head

Security, Safety & Independence

Elevated seats must be secure to be safe. The patented Snap 2 Secure installation system prevents your Clean Shield seat from wiggling or loosening . . . ever!
Easy to install, wrench included.

Why is it important?


A secure seat is a must for anyone with sitting or standing issues
A secure seat provides independence... the user can trust their seat
Loose seats are more difficult to sit on and stand up from, increasing risk of a fall.

CLEAN SHIELD™ elevated toilet seat


PREVENT FALLS: Seat fully secured, will not wiggle or loosen
APPEARANCE: Great looking residential seat that is elevated 3”
SIZE: Available in Round & Elongated
WEIGHT LIMIT: None… tested to 1,000 pounds
SUPPORT ARMS: Available separately… add to seat now or later



Measure your toilet from the center of the mounting hole to the front of the toilet bowl.

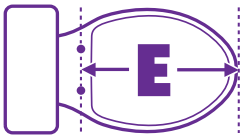

18.9”

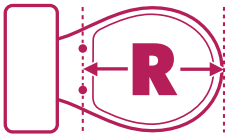

16.5”



CLEAN SHIELD™ Seat

(ROUND) 7YR05300TSS 000

SPECIFICATIONS

- Style: Closed Ring, with Cover
- Size/Shape: Round
- Material: Polypropylene
- Hinges: Hybrid (Stainless steel & polypropylene)
- Hardware Kit: Snap 2 Secure™
- Bumpers: Ring: 4 Super Grip™
- Suggested Uses: Residential & Commercial

MADE IN


Proudly Made in the USA with 75% or More US Content. White Color

PACKAGING & WEIGHTS

- Individual Seat Weight: 6.97 lbs

Individual Packaged Unit:

- Packaged Seat Weight: 9.22 lbs
- Packaging Dimensions: 20.5 x 18.75 x 6.75

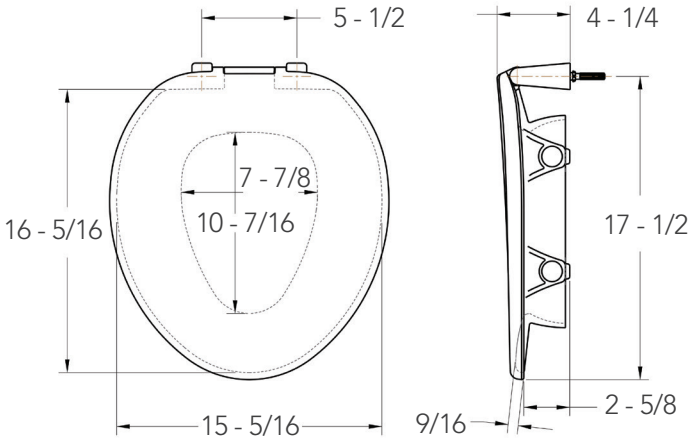

CLEAN SHIELD™ Seat

(ELONGATED) 7YE05300TSS 000

SPECIFICATIONS

- Style: Closed Ring, with Cover
- Size/Shape: Elongated
- Material: Polypropylene
- Hinges: Hybrid (Stainless steel & polypropylene)
- Hardware Kit: Snap 2 Secure™
- Bumpers: Ring: 4 Super Grip™
- Suggested Uses: Residential & Commercial

MADE IN


Proudly Made in the USA with 75% or More US Content. White Color

PACKAGING & WEIGHTS

- Individual Seat Weight: 6.65 lbs

Individual Packaged Unit:

- Packaged Seat Weight: 9.37 lbs
- Packaging Dimensions: 16.37 x 24.5 x 4.18

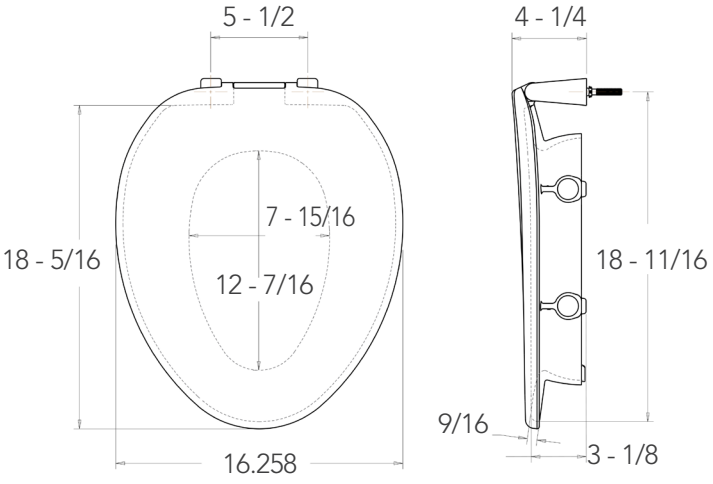

CLEANING: Use warm, soapy water to clean your seat. DO NOT use abrasive cleaners!

SUPPORT ARMS for CLEAN SHIELD™ Seats (7YA05303 GRY)



SPECIFICATIONS

- 21" width between Support Arms
- 1" Support Arm tube diameter
- 1.5" Support Arm handle diameter

MATERIALS


Tube: Heavy Guage Aluminum
Powder coated for protection & appearance
Handles: Polypropylene


WEIGHT CAPACITY

350 pounds each Support Arm

PACKAGING & WEIGHTS

- Individual Seat Weight: 1.8
  Individual Packaged Unit Carton
- Packaged Seat Weight: 2.04 lbs
- Packaging Dimensions: 20.0 x 18.0 x 5.0 lbs

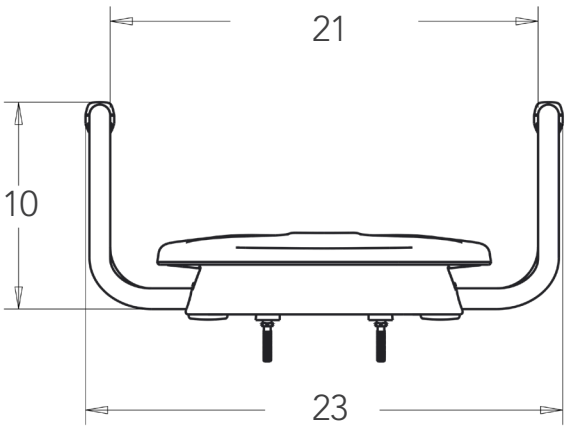

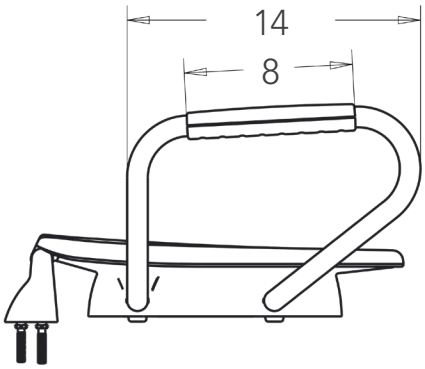

CLEANING: Use warm, soapy water to clean your seat. DO NOT use abrasive cleaners!

SUPPORT ARMS



Support Arms

The Support Arms for the Clean Shield seat are sold separately, so you can add them when you need additional assistance. They are weight bearing, meaning they will support your entire weight. They are exclusively designed for the Clean Shield seat.

How are our arms different?

- Support Arms… full weight bearing, tested to 350 pounds on each side!
- Ergonomic Positioning… proper hand placement helps you use your strength to sit and stand
- Comfortable… 21 inch width between the arms! This width clears the tank, allowing the ring to be placed in the upright position
- Security… Support Arm strength combined with Snap 2 Secure allows the user to push unevenly, or one side only, without any wiggling or loosening of the seat
- Easy to install on your Clean Shield seat





Hinges:

- Clean Shield installs like a regular toilet seat
- Hinges allow ring and cover to be raised for easy cleaning and standing use
- Self-sustaining hinge feature prevents the ring from slamming
- Stainless steel bolts resistant to corrosion

Ring:

- Oversized ring provides additional comfort
- Ring design eliminates pressure points
- Weight capacity: 1,000 lbs.







Seat Cleaning Instructions:

- Toilet seats should be cleaned with mild soap and water… antibacterial hand soap or dish soap can be used
- DO NOT use abrasive cleaners! Abrasive cleaners will scuff and/or scratch plastic toilet seats
- Chemical cleaners are not recommended, but if used, thoroughly wipe down the seat and hinges to remove all cleaner residue

CLEAN SHIELD™ Product Features



NO EMBARRASSING RISER

The Clean Shield elevated toilet seat has a high-end, residential appearance that blends into your bathroom. The 3 inch elevation is virtually unnoticed. The cover provides proper closure when not in use and the smooth, glossy surface matches the ceramic.

CLEAN SHIELD

The innovative "reverse funnel" shield design keeps everything in the bowl. The design of the shield "avoids" contact with waste and the discrete styling of the seat "hides" the shield. When needed, the glossy surface is easy to clean.





KEEPING IT CLEAN

Special elevated hinges on your Clean Shield seat make cleaning very easy. The hinges allow the cover and ring to be placed in the upright position, providing full access to the seat and bowl. The seats are smooth and glossy… easy to clean.

new

4.5” Raised Toilet Seat with Dual Lock Security™

(7YA04505T 000)



A raised seat with support arms safely secured to the toilet with a patented installation system.
Great for:

- Temporary use
- When 3” is not tall enough
- When you need the best “Budget” seat



Features and Benefits include:

Dual Lock Security™

- Two independent locking systems
- Snap 2 Secure attaches the mounting base to the toilet
- Two locking keys attach the seat to the mounting base


Weight Capacity: 300 lbs.

Support Arms

- Included with the seat
- Tested to over 200 lbs. per side

Fits both Elongated and Round Toilets


- Adjustable Mounting Base allows the seat to fit both sizes
- Ring hole is elongated size for additional comfort and personal cleaning access

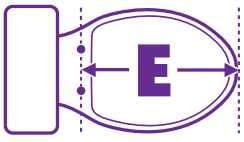

18.9”

Easy to Clean

- Seat has an easy to clean glossy surface
- Unlock to remove for thorough cleaning

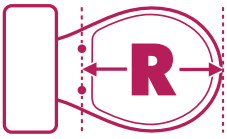

16.5”

new

Personal Wash Bidet™

(7YA06303TWA 000)



Our newest accessory is a bidet attachment specifically designed for the Clean Shield seat. The water stream is calculated to the seating height of the user. The Comfort Stream feature ensures that the water pressure will be "comfortable" for the user.

The bidet installs between the seat hinges and the ceramic using our Snap 2 Secure installation system. The bidet is non-electric, and the water temperature is ambient to the house. There is no easier way to effectively clean oneself.

Personal Wash Bidet™ Features:



- **Comfort Stream™**… Two independent valves control incoming water pressure and water stream strength for widest range of stream options
- **Adjustable T-Valve**… Controls incoming water pressure
- **Handle Control**… Regulates water stream strength when using the bidet
- **Solid Brass Valves**… precision and long life

The Personal Wash Bidet works with all combinations of the Clean Shield seat… Round, Elongated, with or without Support Arms!

Integrity | Trust | Honesty | Knowledge | Respect

WESTECH
HEALTH CARE LTD.

"Our goal to help Canadians of all ages live independently with dignity and respect by offering innovative health care products."



Personal Care Products



Rehab Products

Customer Service



Education & Technical Support



Westech Health Care Ltd.

154 - 5255 McCall Way NE, Calgary, Alberta T2E 7J5
Main: (587) 323-0022 Toll Free: 1 844-323-0022
sales@westechhealth.com
www.westechhealth.com

201120 June 17, 2022 11:57 AM